

Valve XS Atrium Retractor


Updated design providing improved visibility and access



Aesculap Surgical Specialty Products – Cardiothoracic

AESCULAP[®]

Valve XS Atrium Retractor

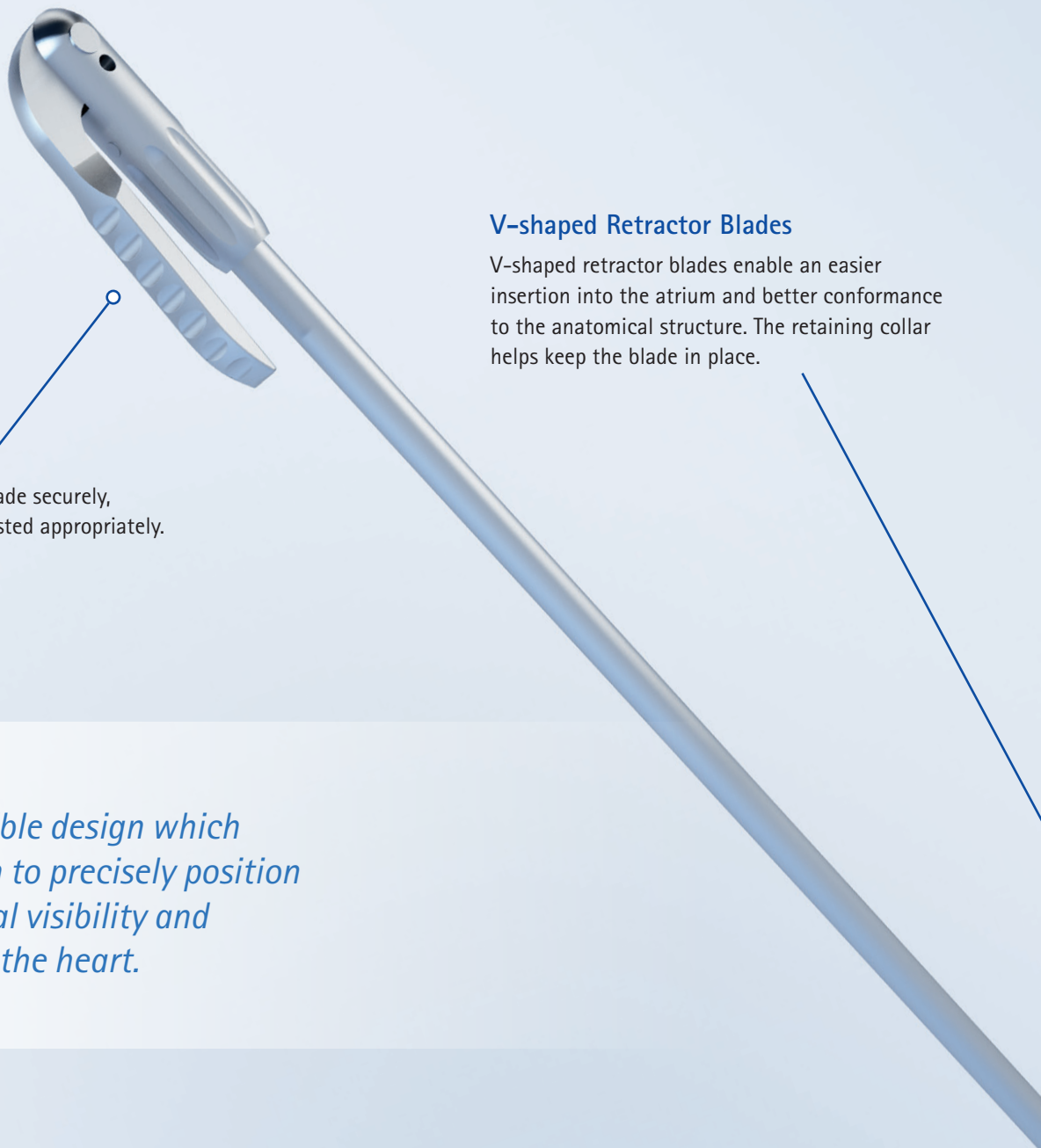

Clamping Lever

Locks the retractor blade securely, after it has been adjusted appropriately.

V-shaped Retractor Blades

V-shaped retractor blades enable an easier insertion into the atrium and better conformance to the anatomical structure. The retaining collar helps keep the blade in place.


Innovative and flexible design which enables the surgeon to precisely position the blade for optimal visibility and increased access to the heart.

Joint Connection

Provides more flexibility in retractor placement, for an excellent mitral valve exposure.

Additional Lateral Blade

Provides an even greater exposure of the mitral valve by retracting more tissue laterally and is easily inserted due to its malleable design.

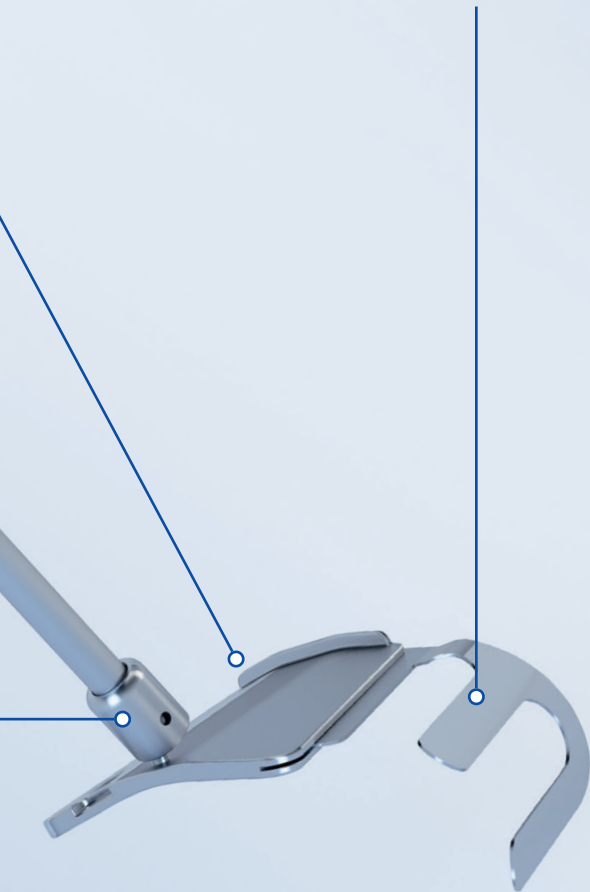

Easy

Easy insertion of the additional lateral blade into the guide groove of the retractor blade.



Quick

Lateral additional blade can be bent to shape, using insertion instrument, after adjustment of the rod and clamping lever.



Better

Excellent mitral valve exposure through additional lateral blade.



Valve XS Atrium Retractor

FC420R

Atrium Retraction Set. Consisting of FC423R - FC430R.

Retractor Blades for Mitral Valve



FC423R

Insertion Instrument
11¾" (300 mm)



FC427R

Rod for Retractor Blades
9½" (242 mm)



Item No.	Height	Width
FC424R	50 mm	35 mm
FC425R	60 mm	35 mm
FC426R	70 mm	35 mm



Item No.	Height	Width
FC428R	50 mm	45 mm
FC429R	60 mm	45 mm
FC430R	70 mm	45 mm

Optional Accessories



2

FC421R

Additional lateral Blade only for Retractor Blades for Mitral Valve (Single use Product)

Blades for Tricuspid Valve



Item No.	Height	Width
FC431R	25 mm	20 mm
FC432R	25 mm	25 mm
FC433R	40 mm	25 mm



Item No.	Height	Width
FC434R	50 mm	25 mm
FC435R	25 mm	40 mm
FC436R	25 mm	50 mm

Mechanical Holding Arm Set, Complete with Atrium Retractor

MG291R Includes Mechanical Holding Arm, OR Table Clamp, Universal Instrument Holder and Atrium Retractor Set (FC420R)



Mechanical Arm



RT090R - OR Table Clamp



RT072R - Universal Instrument Holder

All rights reserved. Technical alterations are possible. This leaflet may be used for no other purposes than offering, buying and selling of our products. No part may be copied or reproduced in any form. In the case of misuse we retain the rights to recall our catalogs and price lists and to take legal actions.

©2015 AESCULAP. ALL RIGHTS RESERVED. PRINTED IN THE USA.
Aesculap is an equal opportunity employer

Aesculap, Inc. | 3773 Corporate Parkway | Center Valley, PA | 18034
Phone 800-282-9000 | Fax 610-791-6886 | www.aesculapusa.com

Aesculap, Inc. - a B. Braun company

DOC1267 Rev A 10/19